

EXPANDING OPTIONS FOR VASCULAR ACCESS

Opus™* Vascular Access Kit

Facilitates serial dilation for peripheral cannulation via the Seldinger technique.



Medtronic
Further, Together

VERSATILITY + CONVENIENCE

Opus™
Vascular Access Kit


**DILATOR
SIZES**
8 Fr to
28 Fr

Stock one
versatile
vascular
access kit to
meet your
needs.

The Opus™ vascular access kit offers a convenient solution for the access you need.

Our kits are thoughtfully composed to meet your preferences and techniques. Each Opus kit includes a complete set of instruments, along with dilators specifically designed to facilitate a smooth insertion — while allowing for gradual dilation of the fascia. A versatile performer, the Opus kit is compatible with single-lumen arterial and venous cannulae — and Crescent™ jugular dual lumen catheters.

- Selection of dilator sizes can accommodate varying patient anatomies.
- Flexible-tip dilators support tracking along the guidewire.†
- Hub helps with handling during insertion and removal.
- Safety scalpel features clear protective shield that locks into place with a one-handed activation.

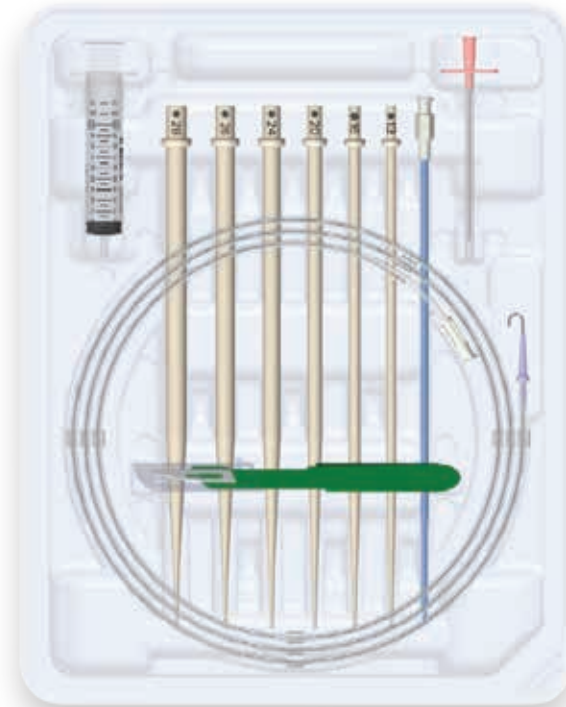

†Data on file at MC3 (200993).



Hubs with
clearly
labeled
sizes

Long tip
gradually
tapers

Ordering Information

CFN	Configuration	Kit Contents	
21030	Box of 5	1 each	Stepped vessel dilator, 8/10 Fr, 13 Fr, 16 Fr, 20 Fr, 24 Fr, 26 Fr, 28 Fr
			Guidewire, J-tip 0.038 in (0.97 mm) x 71 in (180 cm) with depth mark increments
			Scalpel, #11 safety
			Needle, 18 gauge
			Syringe, 10 cc

Experience smooth
insertion and
guidewire tracking
with flexible, long
tapered tip dilators



For more information, contact your local Medtronic Cannula Representative.

For a listing of indications, contraindications, precautions, warnings, and potential adverse events, please refer to the Instructions for Use.

Crescent™ Jugular Dual Lumen Catheter and Opus™ Vascular Access Kit are manufactured by MC3, Inc. and exclusively distributed by Medtronic.

Medtronic

Europe

Medtronic International Trading Sàrl.
Route du Molliau 31
Case postale
CH-1131 Tolochenaz
www.medtronic.eu
Tel: +41 (0)21 802 70 00
Fax: +41 (0)21 802 79 00

United Kingdom/Ireland

Medtronic Limited
Building 9
Croxley Park
Hatters Lane
Watford
Herts WD18 8WW
www.medtronic.co.uk
Tel: +44 (0)1923 212213
Fax: +44 (0)1923 241004

UC201911706EE © Medtronic 2019.
All rights reserved. Printed in Europe.

medtronic.eu